7815 Digitale Personenwaage

de Gebrauchsanleitung en Instructions for Use es Instrucciones de uso fr Mode d'emploi

Inhalt

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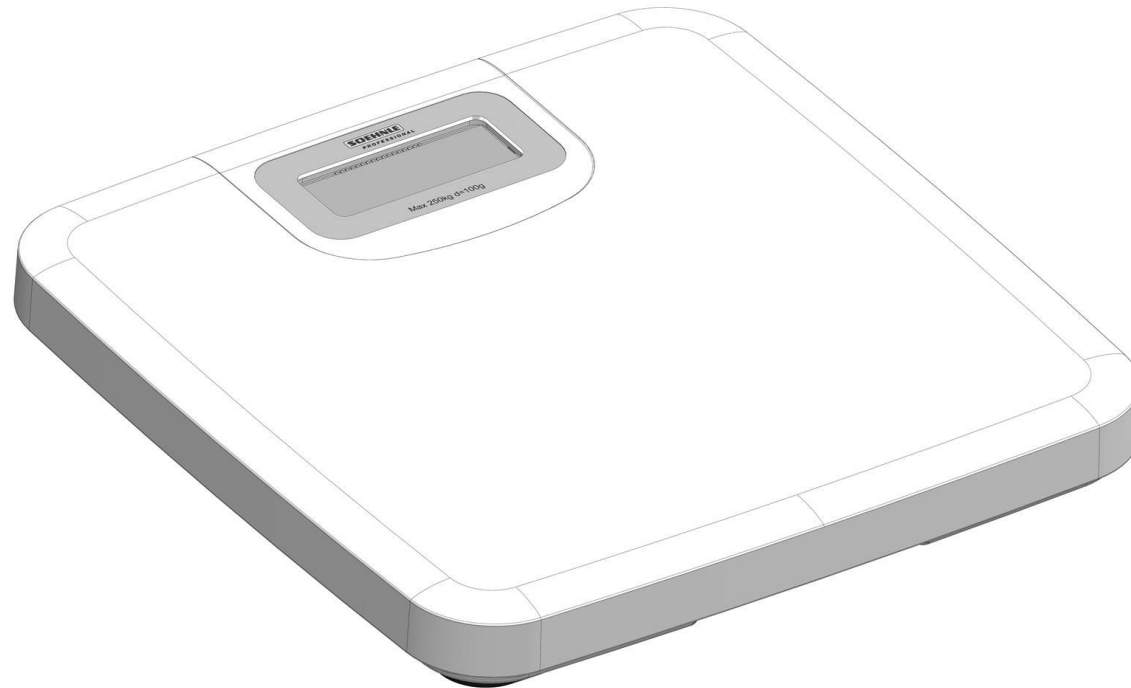
en Instructions for Use

Contents

1. Introduction	1
2. List of parts.....	2
3. Setting up your scale	3
3.1 Setting the kilogram/pound selector switch and installing the batteries.....	3
4.Measuring your weight.....	4
5.Messages on your scale	6
6.Cleaning, maintenance, and disposal.....	8
7.Technical data	8
8. Regulatory.....	10

1. Introduction


The 7815 Digital Personal Scale is a device that you will use to measure your weight.

Read the entire instruction manual before you use your scale. If you have questions about your scale that are not answered by this manual, contact our customer service.


2. List of parts

Ensure that your Soehnle Professional scale package includes the following items:

7815 Digital Personal Scale

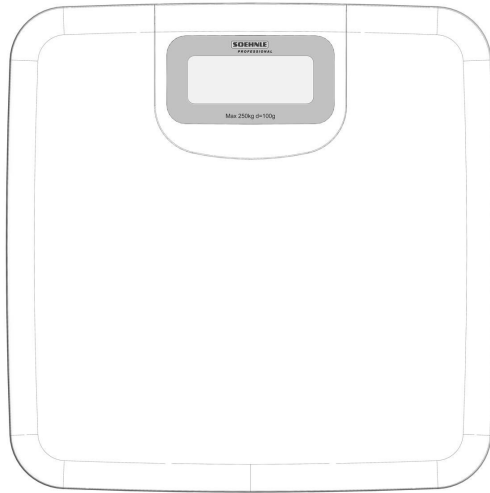Four AA Batteries

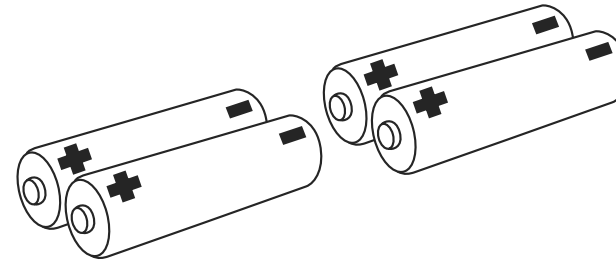Keep the plastic packaging and the batteries away from babies and children – they are a choking hazard.


3. Setting up your scale

Before you use your scale, you need to set the unit of measurement with the kilogram/pound selector switch and install the batteries.

3.1 Setting the kilogram/pound selector switch and installing the batteries

To set the kilogram (kg) / pound (lb) selector switch and install the batteries:

1. Open the battery compartment cover on the bottom of the scale.
2. Inside the battery compartment, set the unit of measurement that is standard for your region by using the **kg / lb** selector switch.

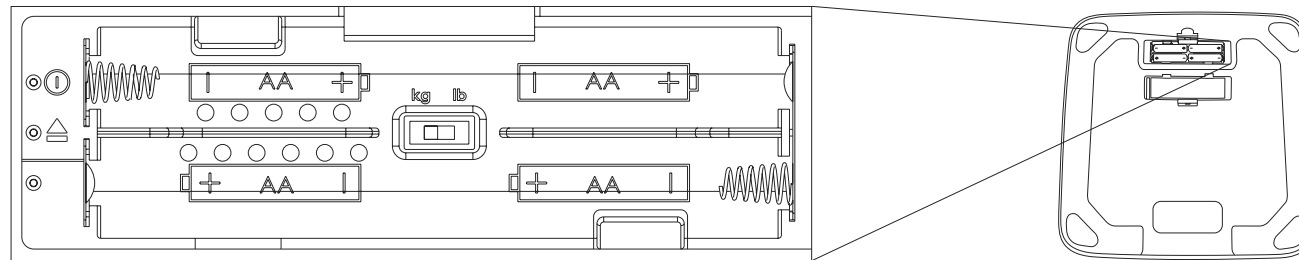3. Insert the batteries as indicated by the diagram inside the battery compartment.
4. Replace the battery compartment cover.

When you replace the batteries, ensure that you replace the complete set of batteries and that you do not use rechargeable batteries, or the battery fluid might leak.

4. Measuring your weight

Regardless of when you measure your weight, the basic procedure is the same.

- Do not step on the scale with wet feet, wet socks, or wet shoes. Ensure that you have stable footing when you step onto the scale.
- Do not let the scale become wet. The scale might be slippery when wet, and liquid might damage the scale.
- Do not support yourself during the weighing process, or your weight measurement will not be accurate.

Before you begin

To ensure accurate weight measurements, observe the following suggestions:

- Try to measure your weight at the same time each day.
- Wear similar clothing each time you measure your weight.

Procedure

To measure your weight:

1. Place the scale on a solid, level, and flat surface.
2. In order to activate the scale, please tap the scale with one foot. The scale is ready when you have heard one beep and 0.0 kg appears on the screen. Step on the scale in order to start the weighing procedure.
3. When your weight is measured, your weight appears on the scale.
4. Step off the scale. The scale displays your weight for approx. 5 seconds and shuts off automatically after 45 seconds.

5. Messages on your scale

When you use your scale, it might display a message. Some messages indicate that a problem occurred.

The following table provides explanations of the messages that you might see on your scale and what you can do to try to resolve problems. If you cannot resolve a problem, contact our customer service.

Message	Number of Beeps	Description
-OL-	5	The weight limit of 250 kg (550 lb) was exceeded.
EC92	5	The calibration data have either been lost or deleted. Please return the scale for inspection.
	0	The batteries are low. Ensure that you have a replacement set of batteries. You can still weigh yourself on the scale.
	5	The batteries must be replaced. Replace the batteries in the scale before you use the scale again.
-UL-	5	Restart the scale unloaded.
H	0	Your weight has been measured.

6. Cleaning, maintenance, and disposal

To clean the scale, use a soft cloth dampened with water. Do not use any chemical or scouring agents. Ensure that liquid cannot penetrate into the scale.

If your scale requires maintenance or disposal, contact your care provider. Do not dispose the scale, accessories, or batteries into household waste.

In Europe, all electrical and electronic waste, batteries, and packaging must be collected separately and disposed of in an environmentally correct manner, according to the following directives: 2012/19/EU.

7. Technical data

The scale has a maximum weight limit of 250 kg (550 lb).

The scale has a minimum weight limit of 8 kg (17 lb).

Use the scale indoors only.

Do not leave objects on the scale. If you leave objects on the scale, the scale might not accurately measure your weight.

If your scale malfunctions, contact our customer service.

Maximum load: 250 kg (550 lb)
Minimum load: 8 kg (17 lb)
Display increments: 0.1 kg (0.2 lb)
Display: LCD, 30 mm character height
Protection class: IP22
Operating temperature: +10°C to +40°C (+50°F to +104°F)
Storage temperature: -30°C to +70°C (-22°F to +158°F)
Batteries: 4 x 1.5 V AA (R6P)

8. Regulatory

Directives

The 7815 Personal Digital Scale scale complies with the relevant regulations of the EU directives and US regulations.

2014/30/EU	EMC
1907/2006/EC	REACH
2011/65/EU	ROHS

Conformity with EU directives

The scale is interference-free in compliance with the applicable EU directive 2014/30/EU. Under extreme electrostatic or electromagnetic influences (for example, when operating a radio unit or mobile phone in the immediate vicinity of the scale) the displayed measurement might be incorrect. When the interference is no longer present, the product can be used again as intended. However, it might be necessary to restart the scale. The scale is a sensitive measuring device. Drafts, vibrations, sudden temperature changes, and sunlight might influence weighing results. The scale complies with protection class IP 22. Avoid exposing the scale to high humidity, vapors, aggressive fluids, and heavy soiling.

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Soehnle Industrial Solutions GmbH